

Padlock Clip™

defect closure system



Full circumferential tissue closure...

The Padlock Clip™ defect closure system facilitates fast and effective full circumferential tissue closure while still maintaining blood flow to promote healing.

The Padlock Clip™ defect closure system features...

- A pre-loaded, self-grasping clip designed to encircle, lift, close, and heal tissue defects
- Quick and easy attachment to the outside of the endoscope
- Open and free instrument channel for optimal endoscope suction
- Intuitive “push of the thumb” deployment

Available in two options:

Padlock Clip™ device for endoscopes between 9.5mm and 11mm

Padlock Clip Pro-Select™ device* for endoscopes between 11.5mm and 14mm



Gastric




Colonic

One clip design for all applications including - hemostasis and closure of ulcers, fistulas, perforations and polyp resection sites < 2cm

*Ability to adjust tissue chamber depth based on scope diameter.

The Padlock Clip™ device is made of flexible, super elastic alloy and lays flat against the tissue, offering low profile radial compression.

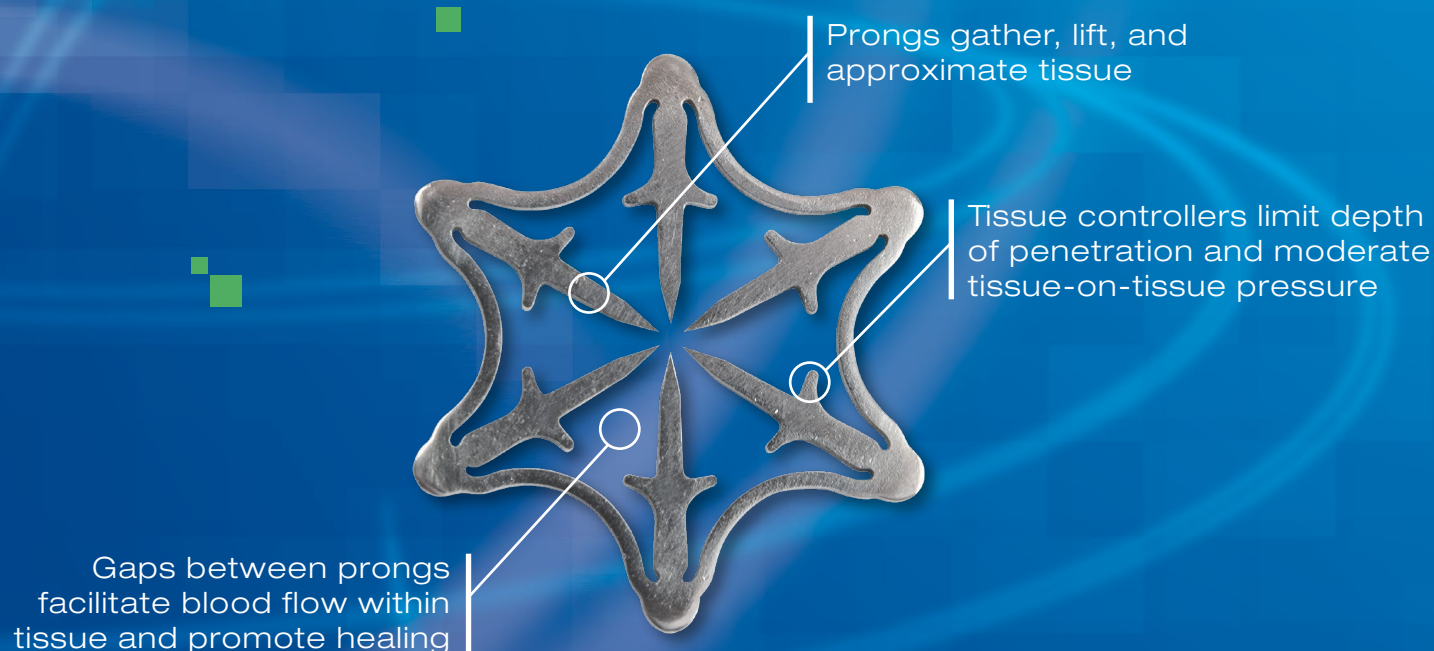

For more information on the Padlock Clip™ defect closure system, contact US Endoscopy today.

product number	description	length (cm)	endoscope tip diameter (mm)	tissue chamber depth (cm)	tissue chamber diameter (cm)	sterile	units/box
C910001	Padlock Clip™ defect closure device	162	9.5-11	1.0	1.1	yes	1
C913131	Padlock Clip Pro-Select™ defect closure device	162	11.5 12.0 12.5 13.0 13.5 14.0	0.4 0.8 1.1 1.3 1.5 1.9	1.1	yes	1



listening... and delivering solutions®

Distributed by:

5976 Heisley Road
Mentor, OH 44060
phone: (1) 440 / 639-4494
fax: (1) 440 / 639-4495
customer service: (1) 800 / 769-8226
email: global@usendoscopy.com

Follow us:    

www.usendoscopy.com

761406 Rev. A