



**DIAGMED**  
HEALTHCARE



## Smartband Ligation System

The new revolutionary device used to endoscopically ligate oesophageal varices at or above the gastroesophageal junction.

### NEXT GENERATION HANDLE DESIGN


Single-hub ligation handle allows for full-functionality with the simplicity of a one-way rotation for band deployment

### PERMACOIL™ LOADING DEVICE

PermaCoil™ technology for reliable device assembly

### COST SAVING REPLACEMENT PACK

When additional bands are required to complete the procedure

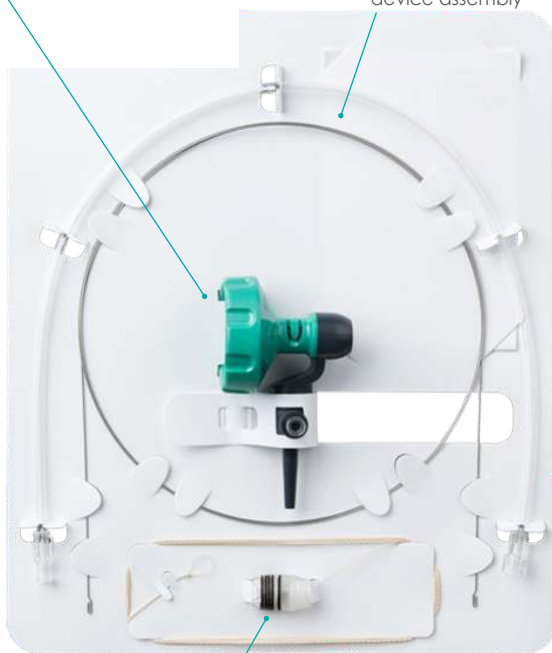

### HIGH RETENTION BANDS

2 X shelf life - twice the effective band life maximising inventory efficiencies

MAKING A DIFFERENCE TO HEALTH

**Features and benefits:**

- High-retention bands contain 50% Latex protein antigens – SmartBand bands contain ½ of the latex protein antigens compared to the current market leading devices. SmartBand has nearly undetectable levels of Latex Protein Antigens
- 100% band retention power
- Specialised UV protected packaging maintains brand performance for upto 24 months shelf life
- Larger handle stem inner diameter, 9.8Fr, ensures maximum suction and flushing capabilities
- Improve endoscopic view during the procedure by using the dedicated flush port
- Revolutionary, single-hub ligation handle allows for full-functionality with the simplicity of a one-way rotation for band deployment
- Enhanced ergonomic control with recessed finger grooves
- The one-piece loading device utilises Perma-Coil technology for reliable device assembly
- Either end of the loading device can be used during setup
- Use a SmartBand Ligation Pack when additional bands are required to complete the procedure



HIGH RETENTION BANDS



FLUSHING PORT

**SmartBand Ligation Kit**

Components: Ligation Handle, Loading Device, Connector Tube and Deployment Cord with Barrel of Bands

Product Code	Barrel Diameter (mm)	Number of Bands	Band Material	Shelf Life (months)
SLK-6	8.6 - 11.6	6	Latex	24

This device is supplied non-sterile and is disposable, single-use only

**SmartBand Ligation Kit**

Components: Deployment Cord with Barrel of Bands

Product Code	Barrel Diameter (mm)	Number of Bands	Band Material	Shelf Life (months)
SLP-6	8.6 - 11.6	6	Latex	24

This device is supplied non-sterile and is disposable, single-use only

