


DuoSwift[®] combination squeegee brush



Combining the cleaning benefits of a bristled brush and squeegee into one device...

The DuoSwift® combination squeegee brush features a tapered nylon brush on the leading end and a series of squeegees on the trailing end to effectively remove debris from flexible endoscopes while adhering to industry guidelines.

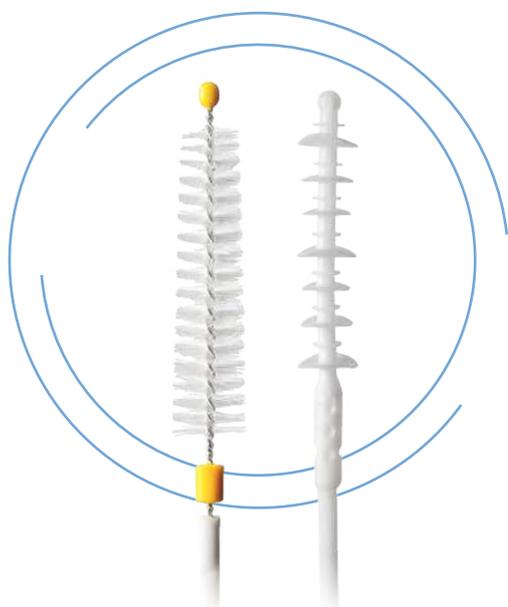

The only combination squeegee brush available

The DuoSwift combination squeegee brush features a...

- » bristled brush to help loosen and remove bioburden with friction
- » squeegee to pull liquid and remaining residue through the endoscope channel
- » distinct design to help maximize cleaning touches within compatible endoscope channels

For more information on the DuoSwift combination squeegee brush, contact STERIS Endoscopy today.

product number	description	sheath diameter (mm)	length (cm)	bristle diameter (mm)	squeegee diameter (mm)	channel size (mm)	units/box
BX00711618	DuoSwift combination squeegee brush	1.7	230	5-7	2.5-5.2	2.8-7.0	50
BX00711619	DuoSwift combination squeegee brush with valve/control head cleaning brush	1.7	230	5-7 10 (valve/control head)	2.5-5.2 (squeegee)	2.8-7.0	50

Distributed by:



© 2020 STERIS Endoscopy or its affiliates. All rights reserved. Unless otherwise indicated, all marks denoted with ® or ™ are registered with the U.S. Patent and Trademark Office, or trademarks owned by STERIS or its affiliates.



STERIS
5960 Heisley Road
Mentor, OH 44060-1834 – USA
Customer Service: 800-548-4873
Fax: 440-639-4495
www.steris.com

Follow us:    