




Snares

As clinical indications for the use of snares continue to grow, so too has our range of innovative, unique and clinically supportive snares.

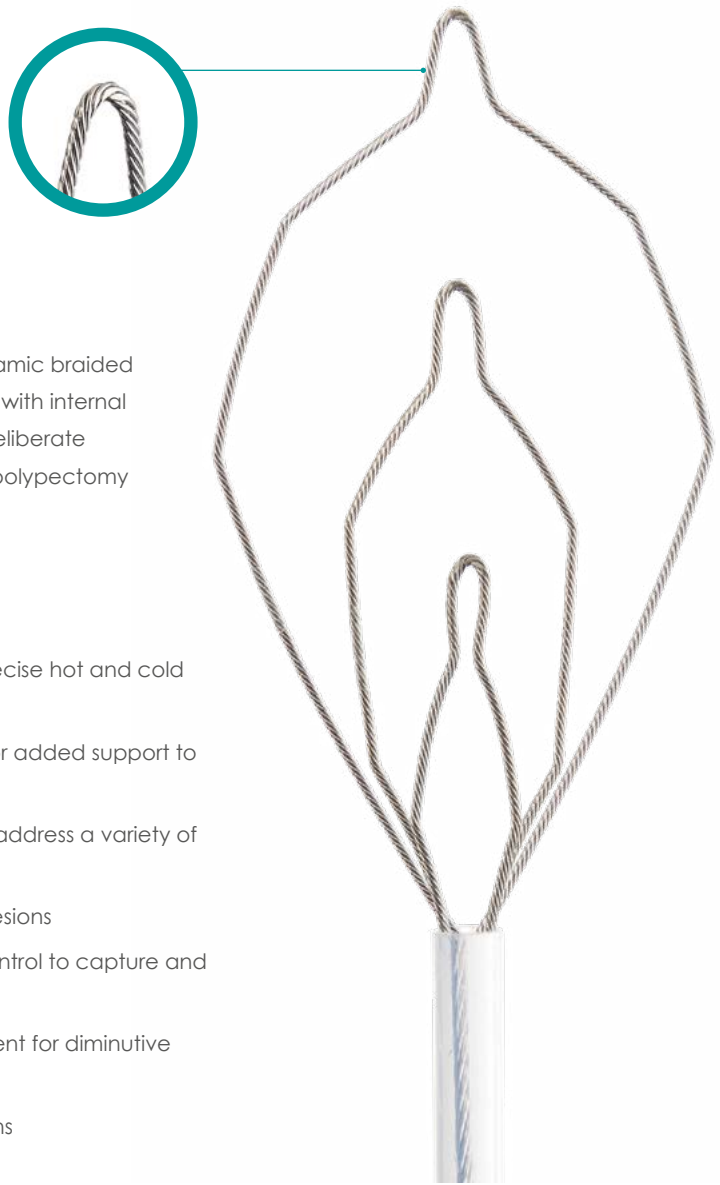

CoinTip™ Snare

The CoinTip™ snare is a hot and cold snare featuring a dynamic braided wire configuration with a coined distal tip, double catheter with internal strain relief, and three distinct shapes. The snare provides deliberate and thorough cutting capabilities specifically for complex polypectomy procedures like endoscopic mucosal resections (EMR).

Features:

- A **flattened, thin wire** at the tip of the snare loop for precise hot and cold polypectomy cutting capabilities
- **Double catheter** that includes an internal strain relief for added support to reduce the chance of buckling during resection
- **Three distinct shapes** and multiple size options to help address a variety of procedural needs
 - » *20MM DELTOID SNARE* - Encompasses majority of lesions
 - » *10MM HEXAGONAL SNARE* - Offers stiffness and control to capture and resect polyps
 - » *6MM DIAMOND SNARE* - Provides specific placement for diminutive polyps and cleaning up resection site margins
- Maintained shape and integrity after multiple resections
- Dual indication for use with or without electrocautery

“The innovative CoinTip Snare supports hot and cold resection with a unique design”



STAGE 1 -
DELTOID 20MM



STAGE 2 -
HEXAGONAL 10MM



STAGE 3 -
DIAMOND 6MM

Benefits of CoinTip™ Snare:


- Provides differentiated solution for complex polypectomies
- Allows en bloc or piecemeal resection abilities with its 20mm size
- Leverages the benefits of different wire diameters to support precise cutting
- Allows technique flexibility with dual indication for both hot and cold use

Ordering Information

| Product Code | Description | Snare Wire Diameter (mm) | Approximate Snare Size (Width x Length) mm | Sheath Diameter (mm) | Sheath Length (cm) | Box Size |
|--------------|----------------|------------------------------------------|--------------------------------------------------------------|----------------------|--------------------|----------|
| 00711110 | CoinTip™ Snare | 0.41 mm (0.016") Tip 0.25 mm (0.010") | Deltoid - 20 x 32 Hexagonal - 10 x 27 Diamond - 6 x 19 | 2.5 | 230 | 10 |

DOUBLE CATHETER

Catheter includes an internal strain relief for added support during resection to reduce the chance of buckling during resection

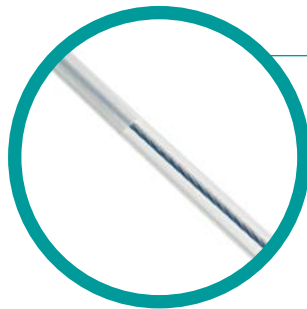

Exacto™ Cold Snare

The unique and exclusive snare design offers more control for placement around diminutive polyps with a precise, clean and safe cut. It is the only custom-made device in the world that is specifically suited for Cold Snaring Techniques.

Features:

- Enables a clean cut, with a snare wire that is 33% thinner than traditional braided snares which makes it better suited to the cold snaring process
- A unique design that is smaller than any available mini snare, providing more control for snare placement and polyp resection
- May help reduce polyp “fly away” from the resection site, due to the unique distal shape



WIRE
33% thinner than
traditional braided wires

“Clinically proven to offer
higher resection rate
than any other snare”



Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|--------------|--------------------|-----------------|---------------|----------------------|-------------|----------|
| 00711064 | Exacto™ Cold Snare | 90 x 190 | Thin, Twisted | 2.4 | 230 | 20 |
| 00711095 | Exacto™ Cold Snare | 90 x 190 | Thin, Twisted | 2.3 | 230 | 50 |

Traxtion™ Hexagonal Snare

Arguably one of the stiffest snares on the market today. The stiff, hexagonal braided wire provides exceptional grip of the mucosa during procedures. The design of the snare incorporates memory points in its nitinol wire to ensure that it consistently opens in the same way multiple times during lengthy procedures.

Features:

- 0.018" diameter wire requiring no changes to electrocautery settings
- Ability to maintain shape and integrity with multiple resections
- Long lasting integrity of the shape may reduce the need for additional snares

“The Traxtion™ Snare is arguably one of the stiffest snares available on the market today”

iSnare®

The iSnare® system is ideal for EMR. This multi-functional device features an alternating injection needle and polypectomy snare which eliminates the need for intraprocedural device exchange, speeding up procedural times.

Features:


- Enables immediate injections in response to bleb absorption/dissipation or bleeding
- Helps maintain scope positioning throughout the procedure which can adjust during instrument exchange, minimising time and effort to relocate the area
- Saves valuable procedure time by lifting, resecting and SPOT® marking polyps without exchanging accessories
- Available in either 23 or 25 gauge needle to suit clinical preference

Ordering Information

| Product Code | Description | Needle | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Sheath Length (cm) | Box Size |
|--------------|---------------------|----------|-----------------|-----------|----------------------|--------------------|----------|
| 00711113 | Traxtion™ Hexagonal | N / A | 25 x 40 | Twisted | 2.3 | 230 | 10 |
| 00711085 | iSnare® Oval | 23G, 5mm | 25 x 40 | Twisted | 3.0 | 230 | 5 |
| 00711086 | iSnare® Oval | 25G, 5mm | 25 x 40 | Twisted | 3.0 | 230 | 5 |
| 00711088 | iSnare® Hexagonal | 25G, 5mm | 25 x 40 | Twisted | 3.0 | 230 | 5 |
| 00711089 | iSnare® Hexagonal | 23G, 5mm | 25 x 40 | Twisted | 3.0 | 230 | 5 |



ISNARE®
Multi-staged, snare types available,
Hexagonal and Oval

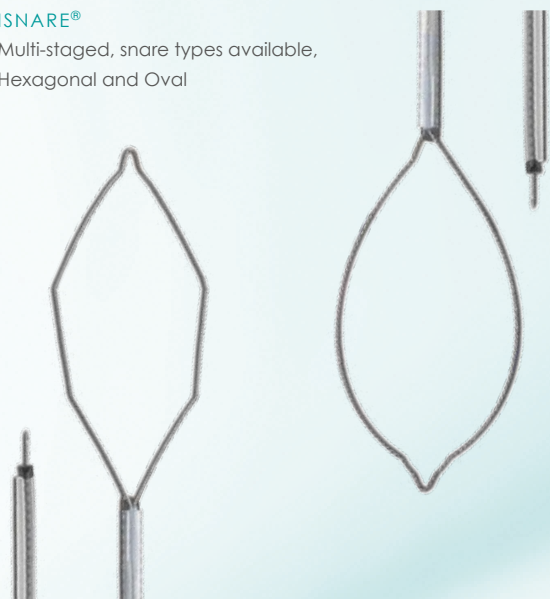

iSnare® Lariat®

The iSnare® Lariat® system is part of a family of unique, multi-functional devices ideal for minimising procedure time and optimising clinical outcomes during submucosal polypectomy procedures, including EMR and SAP.

Features:

- An alternating injection needle and snare operated by a multi-functional handle
- Minimises procedure time by eliminating device exchanges between needle and snare
- Enables immediate snaring to avoid cushion / bleb absorption / dissipation
- Allows immediate injection for haemostasis - if needed after polyp resection
- Features the unique three stage Lariat snare for unrivalled functionality


“Combining a needle and snare into one device to simplify advanced polyp resection...”

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|--------------|-----------------|--------------------------------------------------------|-----------|----------------------|-------------|----------|
| 00711084 | iSnare® Lariat® | Oval: 30 x 55 Hexagonal: 10 x 28 Diamond: 6 x 15 | Twisted | 2.5 | 230 | 5 |



INJECTION
NEEDLE TIP

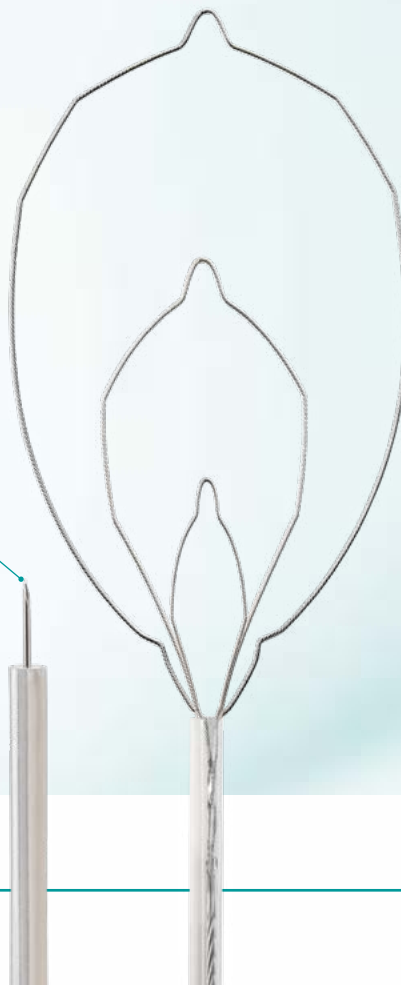



00711084
iSnare® Lariat®

Lariat® Multi Stage Snare

The unique and patented multi-functional design of the Lariat® Multi Stage Snare helps address a variety of clinical needs including piecemeal resection and multiple polyp removal. Three snare shapes in one device can help towards a wide range of clinical requirements when encountered.


Features:

- Three distinct sizes and shapes should eliminate the need to exchange snares
- Snare loop design maintains its shape and integrity after multiple resections
- Stiffer snare wire minimises distal deflection
- Ideal for use when facing varying clinical challenges

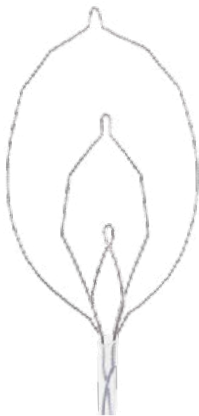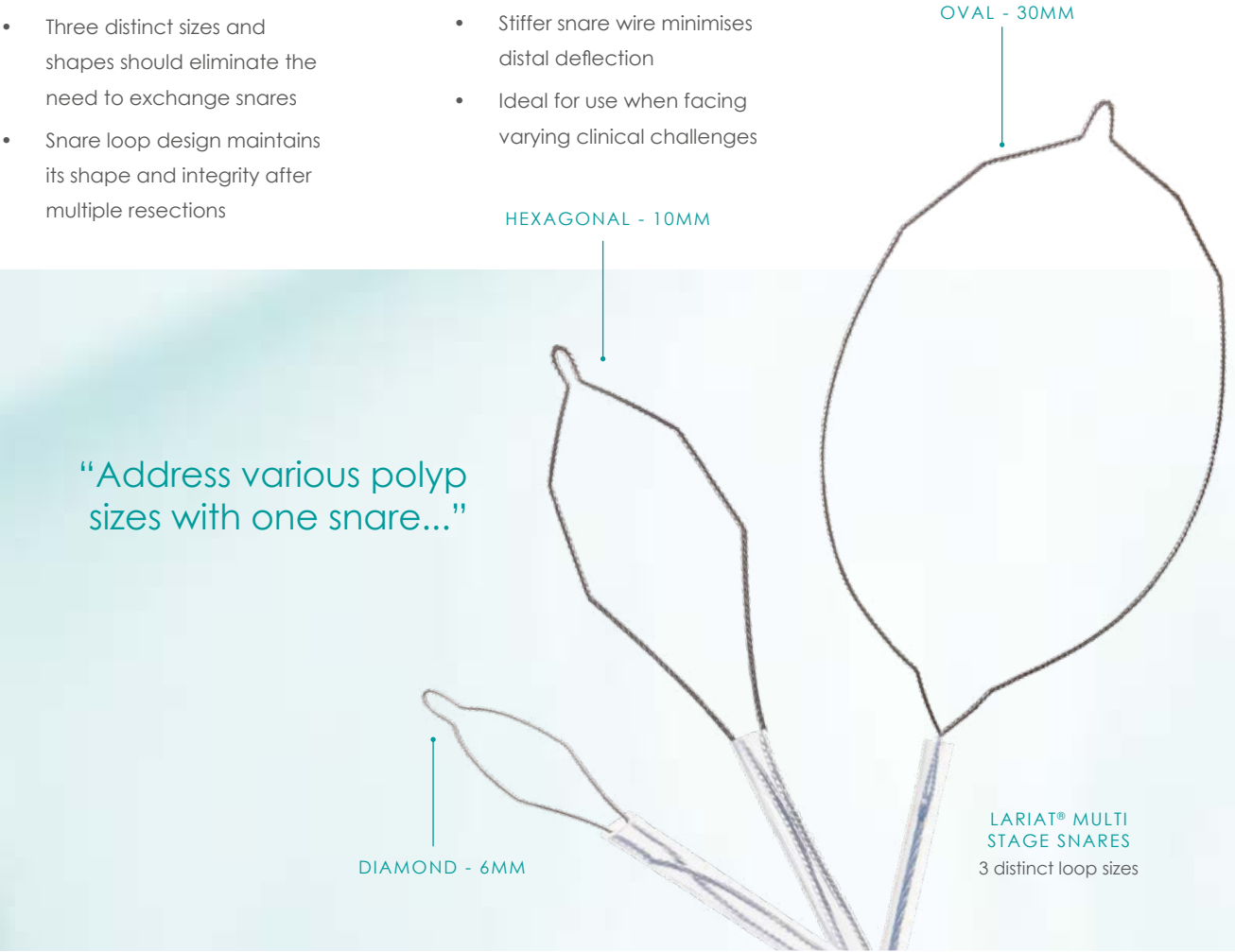

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|--------------|---------------------------|--------------------------------------------------------|-----------|----------------------|-------------|----------|
| 00711119 | Lariat® Multi Stage Snare | Oval: 30 x 55 Hexagonal: 10 x 28 Diamond: 6 x 15 | Twisted | 2.4 | 230 | 10 |
| 00711096 | Lariat® Multi Stage Snare | Oval: 30 x 55 Hexagonal: 10 x 28 Diamond: 6 x 15 | Twisted | 2.4 | 230 | 20 |
| 00711097 | Lariat® Multi Stage Snare | Oval: 30 x 55 Hexagonal: 10 x 28 Diamond: 6 x 15 | Twisted | 2.4 | 230 | 50 |

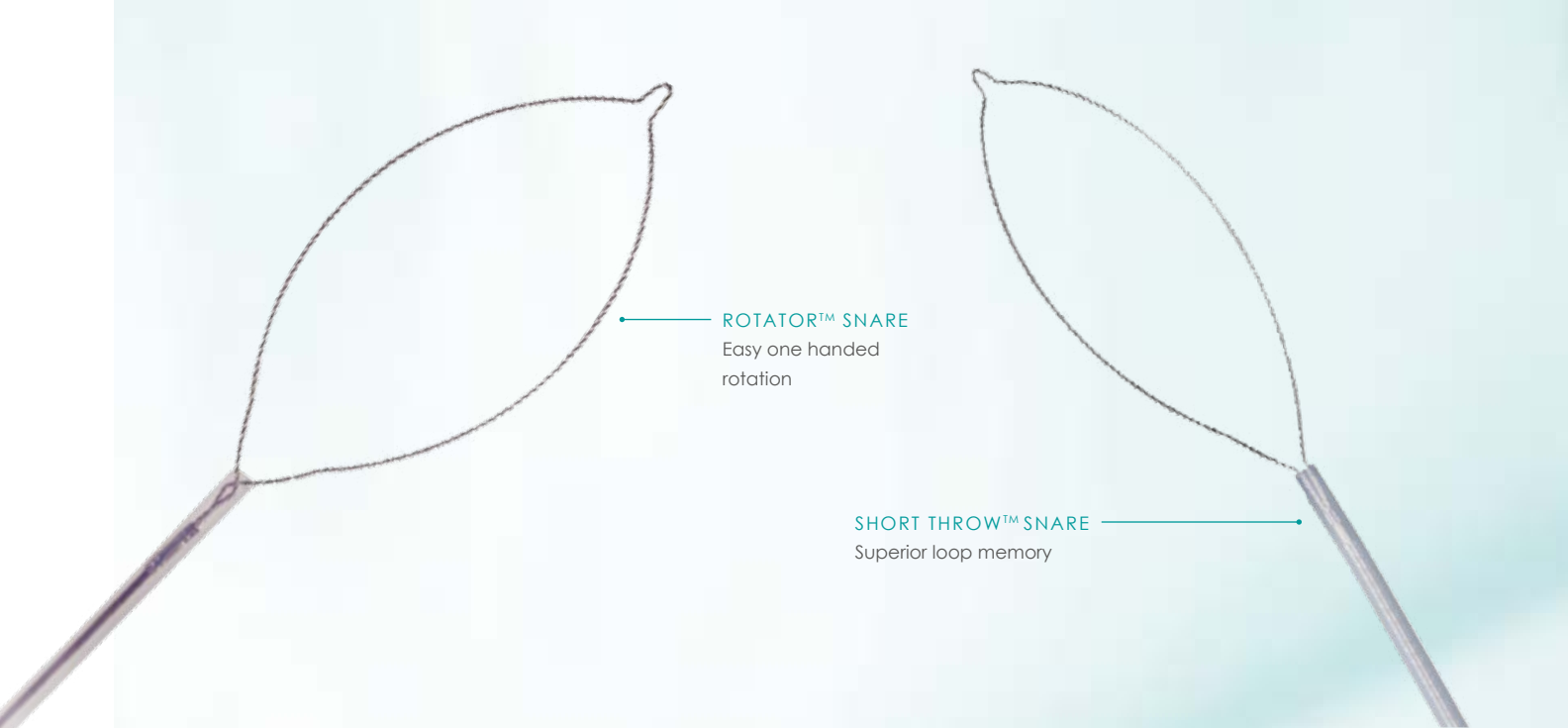

Short Throw™ Snare

The Short Throw™ Snare handle offers a 25% reduction in handle movement, superior loop memory and the comfort of one-handed deployment.

Features:

- The Short Throw™ Snare is designed to work in tortuous conditions; an important feature when performing complex procedures
- Less handle extension is required to fully extend the snare loop
- Suitable as well for those who find non-short throw handles uncomfortable or unusable
- Available in Colonic and Enteroscopic lengths

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Sheath Length (cm) | Box Size |
|--------------|----------------------------------------|-----------------|-----------|----------------------|--------------------|----------|
| 00711018 | Short Throw™ Mini Oval | 12.5 x 30 | Twisted | 2.4 | 230 | 10 |
| 00711019 | Short Throw™ Standard Oval | 25 x 54 | Twisted | 2.4 | 230 | 10 |
| 00711101 | Short Throw™ Standard Oval Enteroscopy | 25 x 54 | Twisted | 2.4 | 350 | 5 |

Rotator™ Snare

This standard or mini oval polypectomy snare easily rotates with one-handed manipulation, which can assist with difficult to reach polyps.

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Sheath Length (cm) | Box Size |
|--------------|-----------------------------------|-----------------|-----------|----------------------|--------------------|----------|
| 00711120 | Rotator™ Standard Mini Oval Snare | 25 x 54 | Twisted | 2.4 | 230 | 10 |
| 00711121 | Rotator™ Mini Oval Snare | 12.5 x 30 | Twisted | 2.4 | 230 | 10 |

Kudo Standard Oval Polypectomy Snare

The Kudo Snare features a controllable loop width opening mechanism which opens up to 45mm to allow the user better control during large flat polyp resection. The dimensions of the snare provide clinical flexibility during EMR procedures.

Features:

- A totally unique snare with variable length and width which uses a stiffer wire
- It is ideal for piecemeal resections where multiple snares can often be needed to complete the procedure

“The Kudo Snare allows the physician better control during large flat polyp resection”



Ordering Information

| Product Code | Description | Snare Opening Size (mm) | Wire Type | Sheath Diameter (mm) | Sheath Length (cm) | Box Size |
|------------------|--------------------------------------|------------------------------|------------------------|----------------------|--------------------|----------|
| 99-05-20-11-23-8 | Kudo Standard Oval Polypectomy Snare | Min: 29 x 57 Max: 45 x 43 | Twisted & Flat Feather | 2.6 | 230 | 2 |

KUDO SNARE

Controllable loop width from 29mm x 57mm through to the widest point of 45mm x 43mm




Polypectomy Snares

ProSnare Polypectomy Snare Range


Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|--------------------|---------------------------------------------------|-----------------|-----------|----------------------|-------------|----------|
| 99-05-20-12-20-2 | ProSnare Oval Braided Wire Polypectomy Snare | 15 x 30 | Twisted | 2.3 | 230 | 10 |
| 99-05-20-12-21-2TW | ProSnare Oval Braided Wire Polypectomy Snare | 20 x 40 | Twisted | 2.3 | 230 | 5 |
| 99-05-20-12-22-2TW | ProSnare Oval Braided Wire Polypectomy Snare | 10 x 20 | Twisted | 2.3 | 230 | 10 |
| 99-05-20-12-23-4TW | ProSnare Hexagonal Braided Wire Polypectomy Snare | 30 x 50 | Twisted | 2.3 | 230 | 5 |

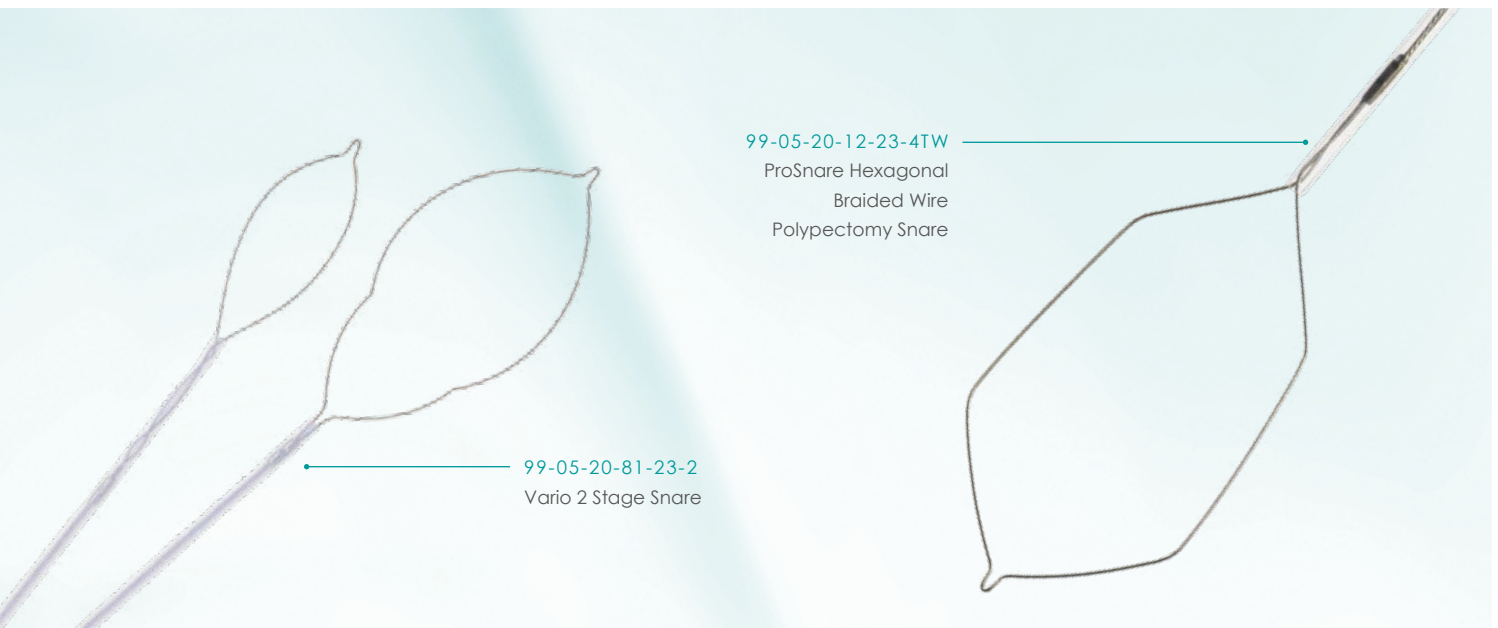

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|------------------|------------------------------------------------|-------------------------------|--------------|----------------------|-------------|----------|
| 99-05-20-81-23-2 | Vario 2 Stage Snare | 15 x 30 30 x 55 | Twisted | 2.3 | 230 | 5 |
| 99-05-21-11-22-5 | Skorpio Polypectomy Snare and Injection Needle | 20 x 40 (needle size: 0.7) | Monofilament | 2.6 | 230 | 1 |
| AG-5071-241523 | Standard Oval Snare | 15 | Twisted | 2.4 | 230 | 1 |
| AG-5073-241523 | Crescent Snare | 15 | Twisted | 2.4 | 230 | 1 |

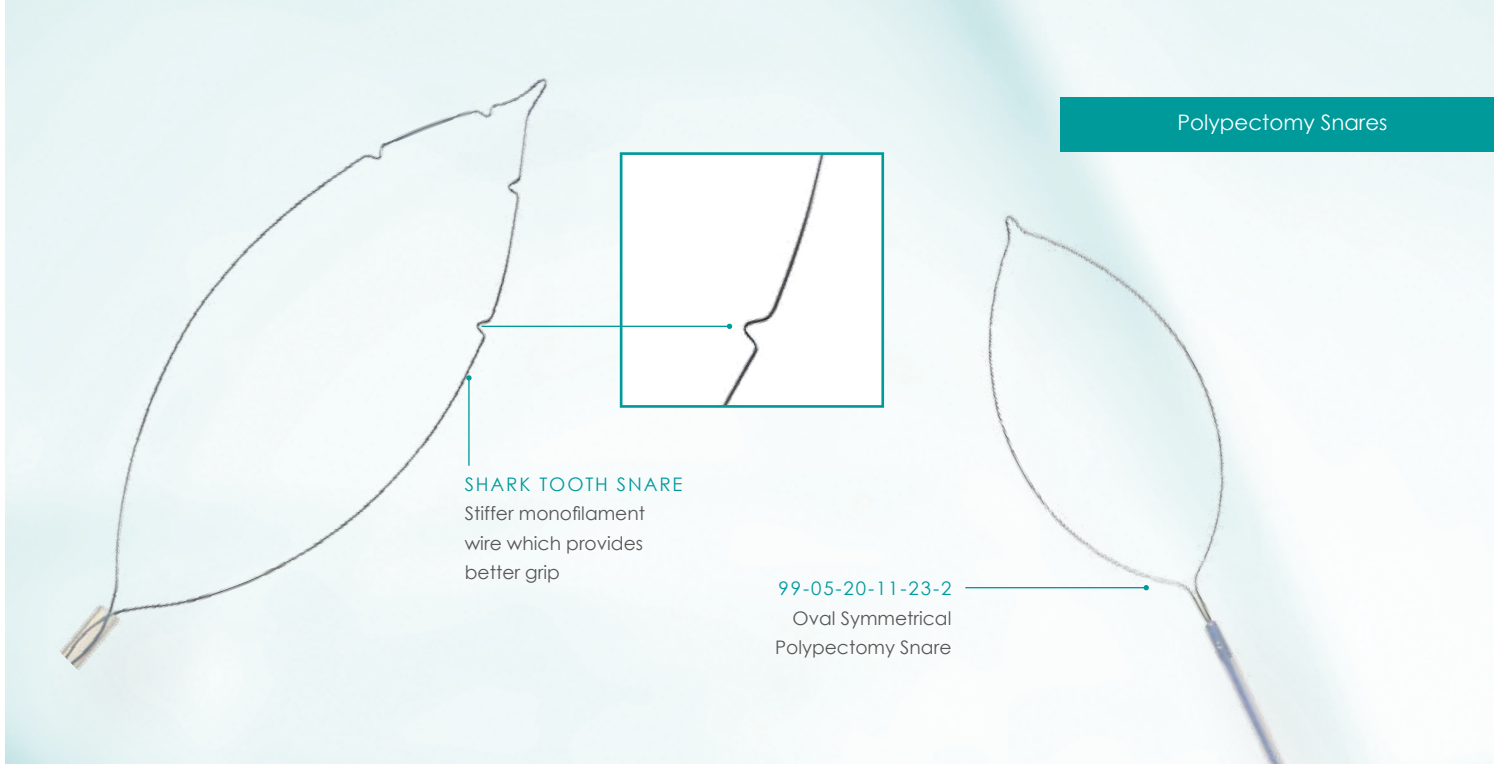

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|------------------|------------------------------------|-----------------|-----------|----------------------|-------------|----------|
| 99-05-20-11-20-2 | Symmetrical Polypectomy Snare | 10 x 25 | Twisted | 2.3 | 230 | 10 |
| 99-05-20-11-21-2 | Symmetrical Polypectomy Snare | 20 x 40 | Twisted | 2.3 | 230 | 10 |
| 99-05-20-11-22-1 | Asymmetrical Polypectomy Snare | 25 | Twisted | 2.3 | 230 | 10 |
| 99-05-20-11-23-2 | Oval Symmetrical Polypectomy Snare | 30 x 55 | Twisted | 2.3 | 230 | 1 |
| 99-05-20-11-24-2 | Symmetrical Polypectomy Snare | 25 x 50 | Twisted | 2.3 | 230 | 10 |

“The asymmetrical snare loop provides the option of safely removing polyps on a fold in the colon”

Shark Tooth Standard Oval Polypectomy Snare

Four unique 'Shark Tooth' shaped barbs located around the distal inner circumference of the wire allow the user to have improved grip and control during resection.

Features:

- The Shark Tooth Snare is manufactured using a stiffer monofilament wire, which provides better snare grip on the mucosa
- The 'shark teeth' grab mucosa and prevent slippage during closure

Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Length (cm) | Box Size |
|-----------------|---------------------------------------------|-----------------|--------------|----------------------|-------------|----------|
| 99-05-20-51-225 | Shark Tooth Standard Oval Polypectomy Snare | 25 x 55 | Monofilament | 2.3 | 230 | 5 |

ezSnare™ Hot Colonic Polypectomy Snare

The ezSnare™ range represents straight forward value for money for non-complex cases and incorporates a short throw handle designed for user comfort that doesn't compromise on effectiveness or performance.

Features:

- Ergonomic handle - with embossed markers to facilitate accurate closure during polypectomy



EZSNARE™
HOT COLONIC
POLYPECTOMY SNARE

Ergonomic Handle
- with embossed
markers to facilitate
accurate closure during
polypectomy



Ordering Information

| Product Code | Description | Snare Size (mm) | Wire Type | Sheath Diameter (mm) | Sheath Length (cm) | Box Size |
|--------------|------------------|-----------------|-----------|----------------------|--------------------|----------|
| SGO-1622S | Small Oval | 10 | Twisted | 2.8 | 220 | 10 |
| SGO-1622M | Medium Oval | 15 | Twisted | 2.8 | 220 | 10 |
| SGO-1622L | Large Oval | 25 | Twisted | 2.8 | 220 | 10 |
| SGO-1622XL | Extra Large Oval | 30 | Twisted | 2.8 | 220 | 10 |
| SGH-1622L | Large Hexagonal | 25 | Twisted | 2.8 | 220 | 10 |



Lumley Close, Thirsk Industrial Park, Thirsk, North Yorkshire, England, YO7 3TD
01845 526660 | info@diagmed.co.uk | diagmed.healthcare

