


Injection Needles



DIAGMED
HEALTHCARE

a STERIS company

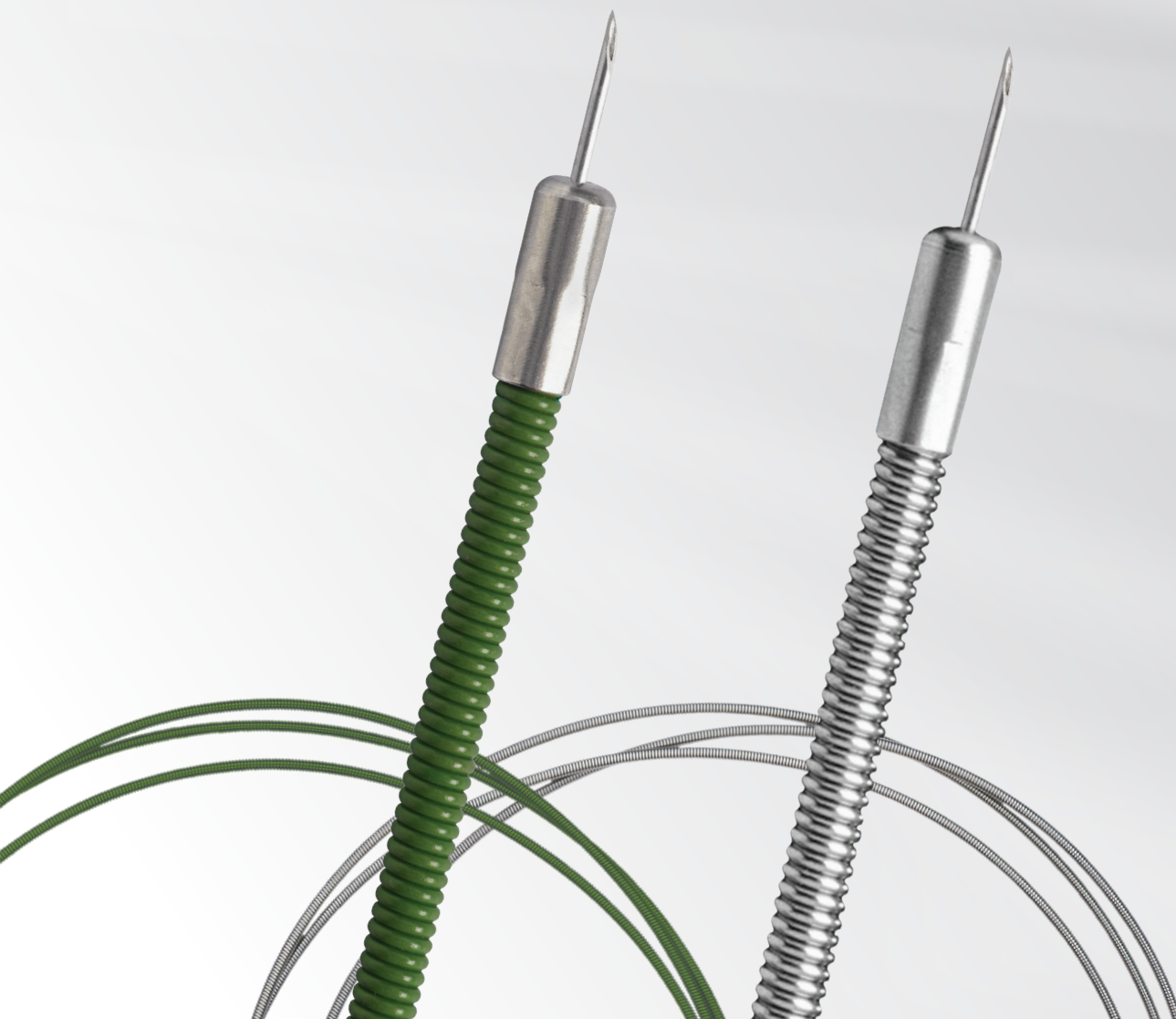

STERIS

Endoscopy

Injection therapy solutions...


STERIS Endoscopy offers the world-renowned Carr-Locke needle and the dependable Articulator™ needle to support your injection needs - from injection assisted polypectomy and enteroscopy to esophageal varices, ulcers, and tattooing.



A green, coiled injection needle with a silver metal hub and a sharp needle tip. It is shown against a light blue circular background. A blue line points from the needle to the right, towards a black endoscope.

Carr-Locke injection needle

Designed with a green Teflon® coating to consistently operate in the most tortuous conditions, including in duodenoscopes

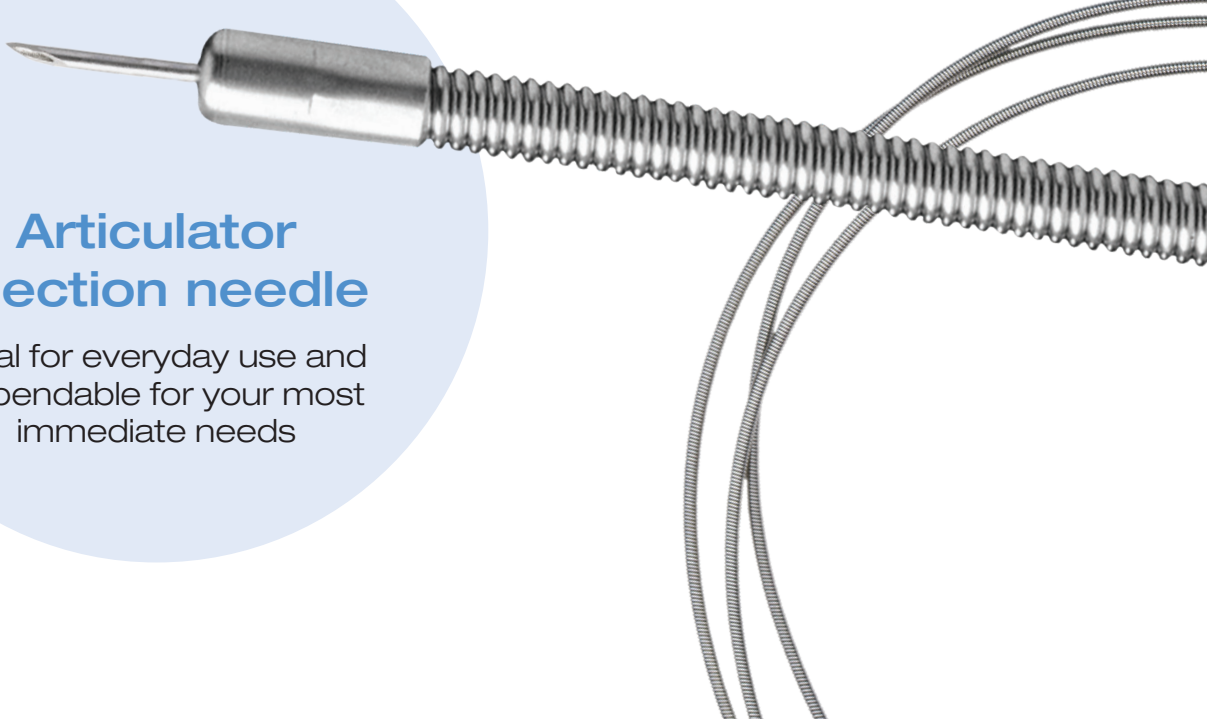


A silver, coiled injection needle with a silver metal hub and a sharp needle tip. It is shown against a light blue circular background. The needle is connected to a black endoscope.

Articulator injection needle

Ideal for everyday use and dependable for your most immediate needs

Features and Benefits

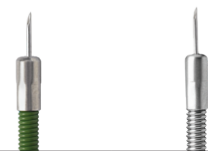

feature	benefit	Carr-Locke needle	Articulator needle
Teflon coated spring sheath and pre-load with stop at distal end	<ul style="list-style-type: none"> Helps to ensure consistent and full length needle projection every time Allows for smooth insertion down endoscope channel 	✓	
Stainless steel spring sheath	<ul style="list-style-type: none"> Minimize sheath kinking especially during challenging procedures 	✓	✓
Ultra smooth distal metal hub	<ul style="list-style-type: none"> Controls depth of needle insertion; reduces the risk of needle perforations Aids in tamponade 	✓	✓
Luer lock, spring loaded handle	<ul style="list-style-type: none"> Automatically retracts needle to reduce risk of inadvertent needle sticks and/or scope damage 	✓	✓
Requires less than 1cc of fluid to prime (most needles require 1.5-2.0cc)	<ul style="list-style-type: none"> Less expensive, particularly when using costly injection agents 	✓	✓

For more information on our injection needles, contact STERIS Endoscopy representative today.

Product Chart

part number	description	sheath diameter (mm)	length (cm)	needle projection (mm)	needle gauge	qty/box
BX00711803	Articulator injection needle	2.5	160	4	25	5
BX00711804	Articulator injection needle	2.5	160	5	25	5
BX00711807	Articulator injection needle	2.5	230	5	25	5
BX00711808*	Articulator injection needle	2.5	350	5	25	5
BX00711810	Articulator injection needle	2.5	230	5	25	10
BX00711811	Carr-Locke injection needle	2.5	230	5	25	5
BX00711812	Carr-Locke injection needle	2.5	230	5	23	5
BX00711813	Carr-Locke injection needle (Pentax compatible)	2.5	230	5	25	5
BX00711814	Carr-Locke injection needle	2.5	230	4	23	5
BX00711822	Carr-Locke injection needle	1.8	230	5	25	5
BX00711823	Carr-Locke injection needle	2.5	230	5	25	10
BX00711824	Carr-Locke injection needle	2.5	230	5	23	10

* BX00711808 does not have the spring-loaded handle

www.steris.com

Distributed by:

PENTAX® is a registered trademark of Hoya Corporation
Teflon® is a registered trademark of E.I. duPont de Nemours and Company as its affiliates

© 2020 STERIS Endoscopy or its affiliates. All rights reserved.
Unless otherwise indicated, all marks denoted with ® or ™ are registered with the U.S. Patent and Trademark Office, or trademarks owned by STERIS or its affiliates.

Follow us:    



Lumley Close, Thirsk
Industrial Park, Thirsk, North
Yorkshire, England, YO7 3TD
01845 526660 |
info@diagmed.co.uk |
diagmed.healthcare